



C-4DL

ICU ELECTRIC TILT BED

Hospital Bed **C-4DL** provides flexible patient care for hospitals that demand safety, security and efficiency in Critical Care and Med Surg. The Profile bed has been independently tested by a globally renowned organization to certify that it complies with the latest stability, dimensional and safety requirements of the IEC60601-2-52 standard.

FEATURES

- ~ Linak® Actuator System, CE, UL safety approved, OpenBus Control Box
- ~ Split Side Rails with Inset Controls for Patient and Nurse
- ~ Retracting Backrest - Reduces Abdominal Squeeze
- ~ Removable HDPE Head & Foot End board, available for CPR
- ~ X-ray box in back panel
- ~ CPR quick release backrest
- ~ One Button chair positioning features
- ~ Modern ABS Deck Panels, Removable and Never-rust (Easy Cleaning & Service)
- ~ Rotating Bumpers at 4 Conner, Soft Contoured Mattress Retainers at Head & Foot
- ~ Integral Base Extension for Tall Patients
- ~ ABS Crank Covers, Dust-proof & easy maintenance
- ~ Battery Back-Up, Under-Bed LED Light
- ~ Imported Tente 6" Twin Central Locking Castors





DIMENSIONS & ADJUSTMENT RANGE

Overall Length	218.5cm
Overall Length Base Extension Extended	238.0cm
Bed Extension Length	19.5cm
Mattress Surface Length - Between Panels	203.0cm
Overall Width	102.0cm
Deck Height Range Twin 150mm Castors	38.5cm to 84.0cm
Trendelenburg Tilt (Electronic)	0 to 14°
Reverse Trendelenburg (Electronic)	0 to 14°
Backrest Angle Adjustment (Electronic)	0 to 65°
Knee Break Angle Adjustment (Electronic)	0 to 35°
Foot Rest Angle Adjustment (Manual)	0 to 28°

RECOMMENDED MATTRESS SIZE DIMENSIONS

Mattress Size	W86cm x L198cm
---------------	----------------

PRODUCT WEIGHT (APPROXIMATE)

With Side Rails Fitted	142Kg
Safe Working Load (Safety Factor = 2xSWL)	250Kg

ELECTRICAL DATA (COUNTRY SPECIFIC)

Power Consumption	230V ac 50Hz
-------------------	--------------



Imported TENTE 6" double castors




DIRECTIONAL LOCK
One direction locking castor. In this position the castor provides better straight line control when the bed is being moved



NEUTRAL LOCK
All functions are released and the bed can be moved freely in any direction



FULLY BRAKED
In this position the bed fully stationary. The profile bed uses four locking castors for complete safety and reliability

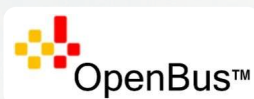

OpenBus Control Box
Brings intelligence to the beds and when combined with our utiliti[®] software solutions, a wide range of asset tracking and service data options can be included.



0 - 65° 0 - 35° 38.5 - 84cm 0 - 14° 0 - 14°

C-4DL

ICU ELECTRIC TILT BED

